
# medisana®

## MG 600

- DE **Massage Gun mit Hot & Cold Funktion**
- GB **Massage Gun with Hot & Cold function**
- NL **Massage Gun met Hot- & Cold-functie**
- FR **Massage Gun avec fonction Hot & Cold**
- ES **Pistola de masaje con función Hot & Cold**
- IT **Massage Gun con funzione Hot & Cold**
- PT **Pistola de massagem com função Hot & Cold**
- GR **Massage Gun με λειτουργία Hot & Cold**
- FI **Massage Gun Hot & Cold -toiminnolla**
- SE **Massagepistol med Hot- och Cold-funktion**
- NO **Massasjepistol med Hot og Cold funksjon**
- DK **Massagepistol med hot & cold-funktion**
- PL **Pistolet do masażu z funkcją Hot & Cold**
- CZ **Masážní pistole s funkcí Hot & Cold**
- SK **Masážna pištoľ s Hot & Cold funkciou**
- SI **Masažna pištola z Hot in Cold funkcijo**
- HR **Pištola za masažu s funkcijom Hot i Cold**
- HU **Massázspisztoly Hot és Cold funkcióval**
- RO **Pistol masaj cu funcție Hot & Cold**
- BG **Massage Gun с функция Hot & Cold**
- EE **"Hot"- ja "Cold"-funktsiooniga massaažipüss**
- LV **Masāžas pistole ar funkcijām "Hot" un "Cold"**
- LT **Masažuoklis „Gun“ su "Hot" ir "Cold" funkcijomis**
- RU **Массажный пистолет с функцией Hot & Cold**

PDF online  
www.medisana.com

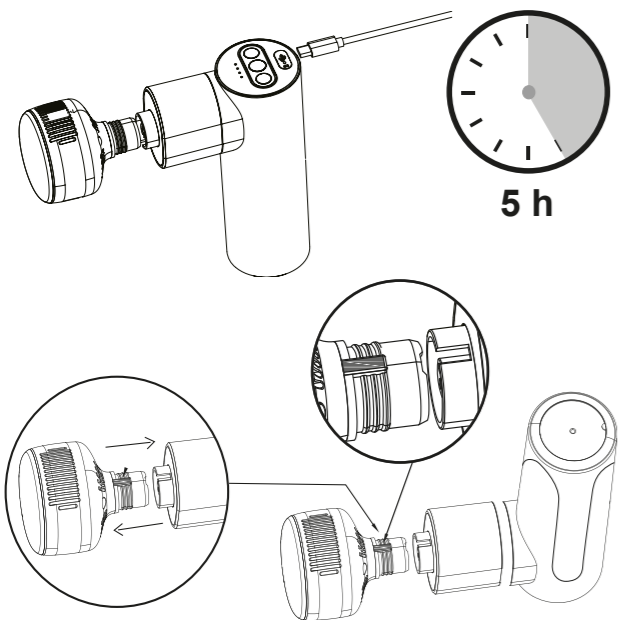5 h

REF 88987

EAN 40 15588 88987 5

medisana GmbH  
Carl-Schurz-Str. 2  
41460 NEUSS, GERMANY

88987 MG600 QSG 08/2022 Ver. 1.0

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**GB Safety Instructions - Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.**

The device is intended for private use only. The appliance must not be used for commercial or medical purposes. · If you have any health concerns, consult your doctor before use. · Use the unit indoors only! Do not use the device in damp rooms (e.g. when bathing or showering). · Use the device only for its intended purpose according to the instructions for use. In case of misuse, the warranty becomes void. · If you feel pain or find the massage uncomfortable, stop the application and consult your doctor. · You must not use the device if you suffer from one or more of the following illnesses or complaints: Circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis or thrombosis. · Do not use the appliance if you are pregnant. Special care is necessary when the appliance is used on, by or near children, sick and helpless persons. · Do not use the device near the eyes or other sensitive parts of the body. · Do not use this device to support or replace medical applications. Chronic conditions and symptoms could worsen. · Do not use the device while you are sleeping or in bed. · Do not use the device before falling asleep. The massage has a stimulating effect. · Do not use the device while driving a vehicle or operating machinery. · This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product in a safe way by a person responsible for their safety. · This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. · Children must not play with the appliance. · Cleaning and user maintenance must not be carried out by children without supervision. · Do not allow the unit to come into contact with pointed or sharp objects. · Do not cover the unit when it is switched on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury. · Do not immerse the unit in water or other liquids. · If liquid does enter the unit, switch it off immediately and disconnect the power supply if it was connected. Wait until the unit is completely dry. · Only connect the unit for charging so that the cable and power connector are freely accessible. · Never leave the unit unattended when it is connected to the mains (for charging) or in operation. · The built-in batteries are not intended for replacement or removal by the end user. If you experience problems with the power supply, contact customer service. · Always switch off the appliance after use by pressing button (Ⓜ). · To disconnect the appliance from the mains, always pull the plug out of the power socket. Do not pull on the cable! · Do not carry, pull or twist the unit by the cable or the power plug. · Make sure that the cable does not pose a trip hazard. Do not bend, pinch or twist the cable. · Do not use the appliance if any damage is visible, if it is not working properly, if it has been dropped or immersed in water. To avoid hazards, send the unit to the service centre for repair. · Do not insert any objects into the openings of the unit. · You may only clean the unit. In case of malfunctions, please do not repair the unit yourself, as this will automatically void all warranty claims. Contact your specialist dealer and have the repairs carried out by the authorised service. · Only use accessories intended for the unit and supplied by the manufacturer. · Excessive use can lead to overheating. After 15 minutes of use, allow the appliance to cool down for **at least** 15 minutes before using it again. The appliance switches off automatically after 15 minutes of use.

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